

Encode® Emergence Healing Transfer

Waarom een Encode®

Natural Tissue Healing

Natural Emergence

- Gradually rising emergence profiles shape and facilitate natural tissue healing.
- Mitigates bone profiling when an implant is placed no deeper than 1-mm below the bone crest.
- Healing abutments are available in a range of emergence profile diameters and in three varying heights: 3, 5, and 7-mm.

Tissue-Level Scan or Impression
Intraoral scan or traditional impression possible when abutment occlusal is clear and above soft-tissue.

Insertion and Removal
Easy abutment insertion and removal as screw remains engaged through screw channel internal threads.

Color and Finish
The pink matte appearance closely matches the gingiva's color for aesthetics results after surgery, and reduces intraoral scan reflection for clear scan images.

MITIGATES THE NEED FOR BONE PROFILING

4-mm (H) Legacy Encode®

3-mm (H)

5-mm (H)

7-mm (H)

Type Encode®

Tapered Screw-Vent® Connection Implants



Platform Diameter	Emergence Profile	Height		
		3	5	7
3.5	3.8	TEEHA333	TEEHA335	TEEHA337
	5.0	TEEHA353	TEEHA355	TEEHA357
	6.5	TEEHA363	TEEHA365	TEEHA367
4.5	4.5	TEEHA443	TEEHA445	TEEHA447
	5.5	TEEHA453	TEEHA455	TEEHA457
	6.5	TEEHA463	TEEHA465	TEEHA467
	7.5	TEEHA473	TEEHA475	TEEHA477
5.7	5.7	TEEHA553	TEEHA555	-
	6.8	TEEHA563	TEEHA565	-
	7.5	TEEHA573	TEEHA575	-

Eztetic Connection Implants

Platform Diameter	Emergence Profile	Height		
		3	5	7
2.9	3.7	CEEHA333	CEEHA335	CEEHA337
	4.5	CEEHA343	CEEHA345	CEEHA347

Encode® Emergence Healing Transfer

Wat zie je

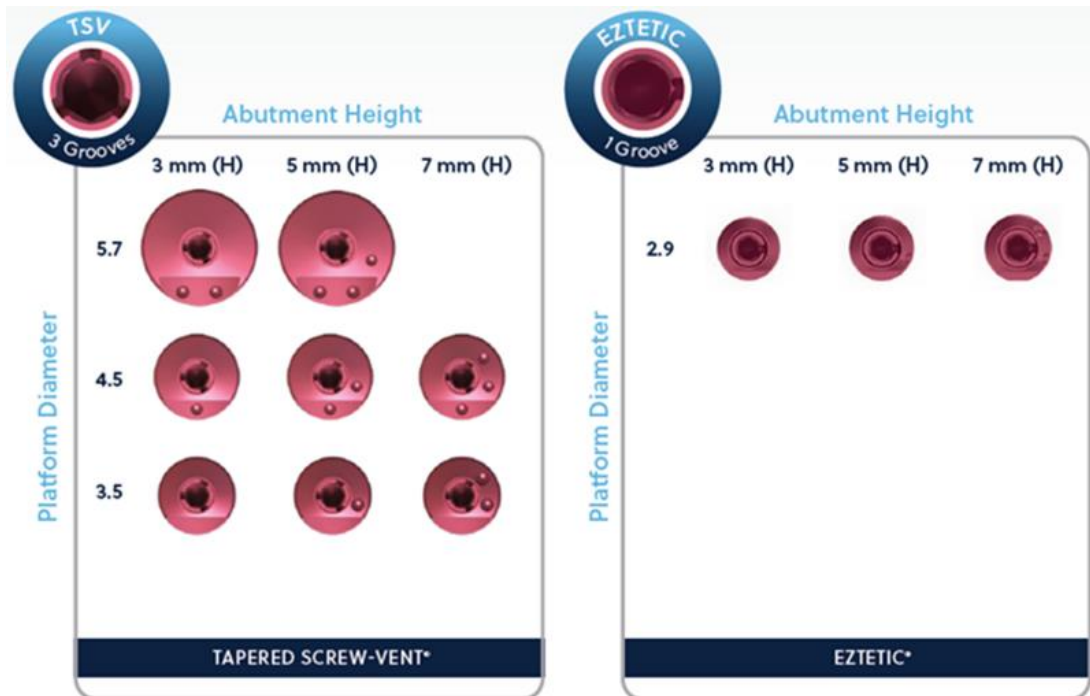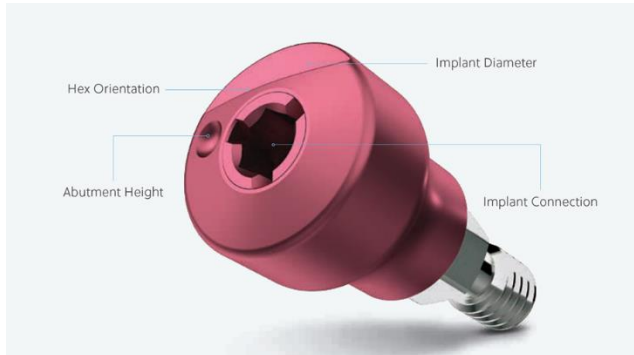

De afdruk

Check altijd of de Encode® correspondeert met het implantaat, type en diameter.

Vermeld naar het laboratorium het type implantaat en type Encode® met diameter, emergence profile en hoogte.

